LigaJect[®]

METRIC

FR

EN

DE

Mode d'emploi	2
Instructions for use	5
Gebrauchsanweisung	8

Your Endo Specialist™

CE

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LigaJect® METRIC

Injection system for intraligamental anesthesia, for use with a 1.8 ml cartridge and metric threaded needle that can be used in intraseptal anesthesia.

EN


INSTRUCTIONS FOR USE

1 Construction of the *LigaJect® METRIC* syringe

The *LigaJect® METRIC* syringe consists of:

- barrel with viewing chamber for anesthetic cartridge
- plastic sleeve
- threaded nozzle for injection needle
- hand piece with dosing lever
- resetting key and dosing plunger.

These two parts are connected by a bayonet catch.

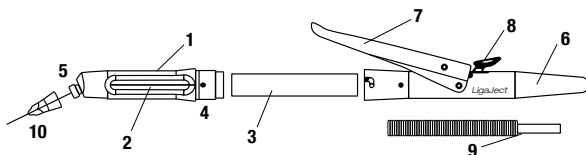- | | | | | | |
|---|--------------------------------|---|--------------------------------------|----|-------------------|
| 1 | Barrel | 4 | Bayonet catch | 7 | Dosing lever |
| 2 | Viewing chamber with cartridge | 5 | Threaded nozzle for injection needle | 8 | Resetting key |
| 3 | Plastic sleeve | 6 | Hand piece | 9 | Dosing plunger |
| | | | | 10 | Disposable needle |

2 Preparation

To return the plunger (9) to the starting point, depress the resetting key (8) and hold the syringe in an upright position (bayonet catch up).

Place the plastic sleeve in the nozzle.

Place anaesthetic cartridge in the barrel and reassemble the two parts of the *LigaJect® METRIC* at the bayonet catch (4).

Depress the dosing lever (7) until a slight resistance is felt. The dosing plunger is now in contact with the anaesthetic cartridge bung and presses the tip of the cartridge into the inner part of the nozzle.

Insert a disposable needle into the nozzle (5) and screw the hub onto the thread. We recommend the use of the cartridge of disposable needles in sizes 0.3 x 12 mm or 0.3 x 16 mm with the corresponding short needle point.

Depress the dosing lever until the anaesthetic solution emerges from the needle point. The syringe is ready for use.

LigaJect® METRIC is sold non-sterile. Prior to first use, carry out pre-disinfection and sterilization. (See paragraph 5)

3 Handling the **LigaJect® METRIC** syringe

Each depression of the dosing lever applies approximately 0.06 ml of anesthetic.

Normally 0.15 ml to 0.20 ml of anesthetic is required for each root, i.e. 3 depressions of the lever.

Multi-rooted teeth need an appropriate dose of anesthetic.

4 Performing an intraligamental anesthesia with **LigaJect® METRIC**

Intraligamental anesthesia provides effective and immediate anesthesia of a single tooth.

The recommended procedure is as follows:

- The needle point is apically inclined towards the central axis of the tooth.
- In contact with the tooth the point of the needle is slowly inserted approximately 1 - 2 mm into the periodontal ligament.
- The needle is correctly placed in the periodontal space if a steady, uninterrupted back-pressure is noticeable.
- This resistance must be maintained, to secure penetration of the solution to the apex.
- Inject slowly and uniformly.

One depression of the dosing lever should take about 7 seconds, i.e. for 0.18 ml per root = 3 depressions of the dosing lever = 21 seconds.

The onset of anesthesia is generally obtained within the 30 seconds following the injection.

Hence the treatment can be started at once.

5 Cleaning and care of the **LigaJect® METRIC**

Please follow these instructions to maintain your **LigaJect® METRIC** syringe in perfect working condition:

- After use, separate the hand-piece and barrel at the bayonet catch. Remove the protective anesthetic cartridge, plastic sleeve and dosing plunger. The **LigaJect® METRIC** may be pre-disinfected or soaked with any product compatible with stainless steels, whilst respecting the conditions of time and concentration indicated by the supplier or manufacturer of the product.
- **Cleaning: the stainless steel LigaJect® METRIC syringe can be cleaned together with all other stainless steel instruments using any kind of common cleaning method at the dentist's.** Immediately after the immersion phase, rinse the parts of the syringe thoroughly, the dry them meticulously.
- **Sterilization:** Place the parts of the **LigaJect® METRIC** in a hermetic sterilization bag. **LigaJect® METRIC** can be sterilized in an autoclave up to 135°C / 2 bars (275°F / 29psi) for 18 min. **The stainless steel LigaJect® METRIC can be sterilized after each use.** Open the sachet just before use on the patient.
- Always ensure that components are completely dry before re-assembly. No need to oil.
- Use only cartridge with dosing plunger tip without cavity.

WARRANTY

Stringent quality control procedures ensure this surgical instrument is supplied in perfect working order. Your **LigaJect® METRIC** is guaranteed against defects in materials or workmanship for **ONE YEAR**

from the date of purchase. This guarantee is void if your **LigaJect® METRIC** has been improperly sterilized, poorly stored or incorrectly handled.

If you return your **LigaJect® METRIC** syringe, make sure you include the warranty card, and please send it to the dental dealer from which it was purchased.

The warranty will only operate if this warranty card is completed.

Therefore, please clip the warranty card onto the accompanying invoice from your dental depot.

Our technical advice, whether provided verbally, in writing or by way of trials, is given in good faith but without warranty, and this also applies where the proprietary rights of third parties are involved. It does not release you from the obligation to test the products supplied by us as to their suitability for the intended processes and uses. The application, use and processing of the products are beyond our control, and therefore entirely your own responsibility.

Should, in spite of this, liability be established for any damage, it will be limited to the value of the goods delivered by us and used by you. We will, of course, provide consistent quality of our products within the scope of our General Conditions of Sale and Delivery.

USA: Caution: Federal law restricts this device to the sale by or on the order of a dentist.

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